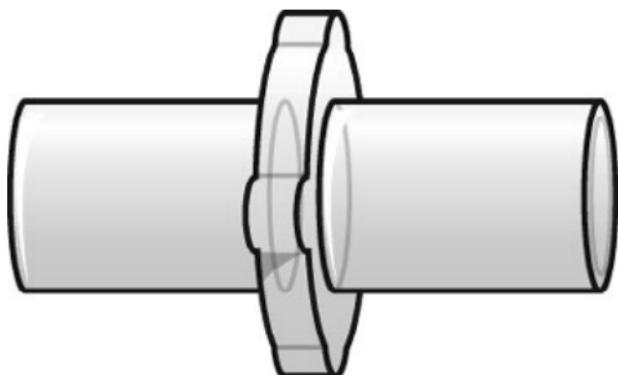
Adaptor 10 m-10 m

**USER INSTRUCTIONS**

**REF 900MR143**

---

<b>English</b> .....	3
<b>Български (Bulgarian)</b> .....	5
<b>Hrvatski (Croatian)</b> .....	7
<b>Česky (Czech)</b> .....	9
<b>Dansk (Danish)</b> .....	11
<b>Nederlands (Dutch)</b> .....	13
<b>Eesti keel (Estonian)</b> .....	15
<b>Suomi (Finnish)</b> .....	17
<b>Français (French)</b> .....	19
<b>Deutsch (German)</b> .....	21
<b>Ελληνικά (Greek)</b> .....	23
<b>Magyar (Hungarian)</b> .....	25
<b>Italiano (Italian)</b> .....	27
<b>Latviešu (Latvian)</b> .....	29
<b>Lietuvių (Lithuanian)</b> .....	31
<b>Norsk (Norwegian)</b> .....	33
<b>Polski (Polish)</b> .....	35
<b>Português (Portuguese)</b> .....	37
<b>Português (Brasil) (Brazilian Portuguese)</b> .....	39
<b>Română (Romanian)</b> .....	41
<b>Slovenčina (Slovak)</b> .....	43
<b>Slovenčina (Slovenian)</b> .....	45
<b>Español (Spanish)</b> .....	47
<b>Svenska (Swedish)</b> .....	49
<b>Türkçe (Turkish)</b> .....	51
<b>繁體中文版 (Traditional Chinese)</b> .....	53
<b>Bahasa Indonesia (Indonesian)</b> .....	55
<b>ไทย (Thai)</b> .....	57
<b>Tiếng Việt (Vietnamese)</b> .....	59
<b>Русский (Russian)</b> .....	61
<b>한국어 (Korean)</b> .....	63
<b>日本語 (Japanese)</b> .....	65

**REF 900MR143****Adaptor 10 m-10 m****Intended Use**

The 900MR143 adaptor is used to connect an inspiratory extension tube to the patient end temperature probe adaptor for neonatal patients in incubators. This accessory is intended for use with Fisher & Paykel Healthcare 900MR780, 900MR781, and 900MR782 neonatal reusable breathing circuits in hospital or clinical environments and should only be prescribed by a physician.

**Product Specifications**

<b>Interface Connections</b>	ISO 5356-1 conical connectors
<b>Maximum Operating Pressure</b>	80 hPa (cmH <sub>2</sub> O) (8 kPa)

**Warning, Cautions and Notes****WARNINGS**

- Breathing circuits and accessories should be prescribed and initially set up by a healthcare professional qualified in respiratory care and used under the supervision of trained medical personnel.
- Do not lubricate connections, tubing or accessories. Failure to comply may result in fire, burns and loss of system performance.
- This product is reprocessable (reusable). It must be reprocessed according to the reprocessing instructions for this product. Incorrect reprocessing may lead to the risk of cross-infection between patients and potential serious harm.
- Visually inspect products before each use on a patient. Discard if faulty or if there is any sign of deterioration such as cracks, tears or damage. These may cause gas leaks and loss of ventilation or respiratory support.
- Check all connections are tightly secured before use.
- Use of non-compatible, non-ISO 5367, or non-ISO 80601-2-74 compliant circuit components may cause unintentional tube disconnection and loss of ventilation or respiratory support.
- Do not modify this product.
- Do not use breathing circuits, chambers, accessories or combinations which are not approved by Fisher & Paykel Healthcare.

**Failure to comply with the above warnings and cautions may impair the performance of the product or compromise safety (including causing potential serious harm).**

- California residents, please be advised of the following, pursuant to Proposition 65: This product contains chemicals known to the State of California to cause cancer, birth defects, and other reproductive harm. For more information, please visit: [www.fphcare.com/prop65](http://www.fphcare.com/prop65)

**Notes**

- Before use on a patient, all components must be cleaned and disinfected. Sterilization is highly recommended if products can withstand the process.
- Ensure the complete circuit functions correctly, with the required ventilator settings, before connecting to a patient.
- This product is designed for the delivery of air, oxygen or a mix of both gases. It is not suitable for the delivery of flammable anesthetic gas mixes, Heliox or gas solutions, suspensions or emulsions that have not been evaluated.
- The responsible organization is accountable for the compatibility of the humidifier and all of the parts and accessories used.
- This product is safe for use on children.

If a serious incident has occurred while using this device, please contact your Fisher & Paykel Healthcare representative and Competent Authority.

Contact your local Fisher & Paykel Healthcare representative for a printed copy of these user instructions.

## Reprocessing Instructions

For reprocessing instructions (UI-623716), please visit [www.fphcare.com/fp-reusables](http://www.fphcare.com/fp-reusables)

### **WARNINGS**

- The product must be disconnected from both the power and gas source prior to cleaning.
- Clean the product prior to use and whenever the product is visibly soiled. Follow the reprocessing instructions (UI-623716). Alternative cleaning methods may deteriorate or damage the product, reducing the product life and compromising therapy.
- Discard the product after 50 reprocessing cycles, or five years from date of manufacture, whichever occurs earlier.

### **Notes**

- It is recommended that reprocessing commences as soon as is reasonably practical following use to prevent gross contaminants from drying on the product.
- Used products must be treated as being contaminated. The user may be exposed to breathing tract fluids. Follow all local, state and federal regulations with respect to environmental protection and correct disposal of these products.

### **Storage**

- As a minimum, reprocessed breathing circuits and accessories are to be stored in a clean plastic bag or container and away from direct sunlight or heat if not being reused immediately.

## Symbol Definitions

		
Product does not contain phthalates or pyrogens	Transportation and storage temperature limits	European Union authorized representative
		
Not made with natural rubber latex	Caution symbol for warnings	Reference number
		
Prescription only	Type BF Applied Part	Date of manufacture
		
Instructions for Use	Do not recycle	Lot number
		
Use-by date	CE Mark	Manufacturer
		
Medical Device		

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