
Nitric Oxide Adapter Kit

USER INSTRUCTIONS

REF 950X08

CE 0123 Single Use

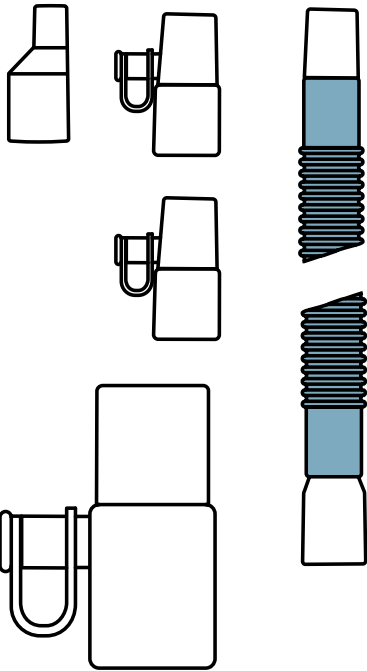English en

F&P 950™ Nitric Oxide Adapter Kit

Indications for use

The adapter set is an accessory to the F&P 950 Respiratory Humidifier. It is intended for use with an F&P 950 breathing set.

Contraindications

There are no contraindications for heated respiratory humidification.

Side effects

There are no known side effects of heated respiratory humidification.

Technical specifications

Compatible with the F&P 950 Respiratory Humidifier. Refer to humidifier, breathing set, and nitric oxide delivery system user instructions for additional warnings, cautions, contraindications, side effects, and system information.

Interface connections: ISO 5356-1 Conical Connectors

Setup

Nitric oxide to be administered and monitored as per hospital protocol.

Symbol definitions

Follow instructions for use - safety.

Consult instructions for use. fphcare.com/950IFU

Single use.

Type BF applied part.

This product is not made with natural rubber latex.

Temperature limitations.

Manufacturer.

Date of manufacture.

Date of expiration.

Warnings, cautions, and notes

WARNINGS

- Appropriate patient monitoring (e.g. oxygen saturation) must be used at all times. Failure to monitor the patient (e.g. in the event of an interruption to gas flow) may result in serious harm or death.
- Single use. Do not reuse this product. Reuse may result in infection of the airway.
- The use (or modification) of breathing circuit/chamber combinations not recommended by Fisher & Paykel Healthcare may result in poor humidification, ventilator malfunction, or serious harm.
- Do not use product near any ignition source. Exposure to oxygen increases the risk of fire that may cause serious injury or death.

Failure to comply with the following warnings may impair performance of the device or compromise safety (including potentially cause serious harm):

- This product is designed for the delivery of air, oxygen, and/or nitric oxide. It is not suitable for the delivery of flammable anesthetic gas mixes or Heliox gas.
- Do not crush, stretch, or milk the tubing.
- Set appropriate ventilator or flow source alarms to monitor therapy delivery.
- Perform a pressure and leak test on the breathing system before connecting to a patient.
- Monitor circuit condensate to prevent occlusion or build-up of fluid. Drain as required.
- Do not use beyond the 14-day maximum duration of use.
- Check all connections are tight before use.
- Do not clean or sterilize this product. Avoid contact with chemicals, cleaning agents, or hand sanitizers.

NOTES

- User may be exposed to breathing tract fluids during disposal. Breathing circuits must be treated as contaminated. Dispose of according to local, state or country specific environmental laws and regulations.
- If a serious incident has occurred while using this device, please inform your local Fisher & Paykel Healthcare representative and, for European Union member countries, the Competent Authority in your country.

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