
F&P Optiflow[™] Nasal Interface with CO₂ Sampling

REF AA030J S/M/L

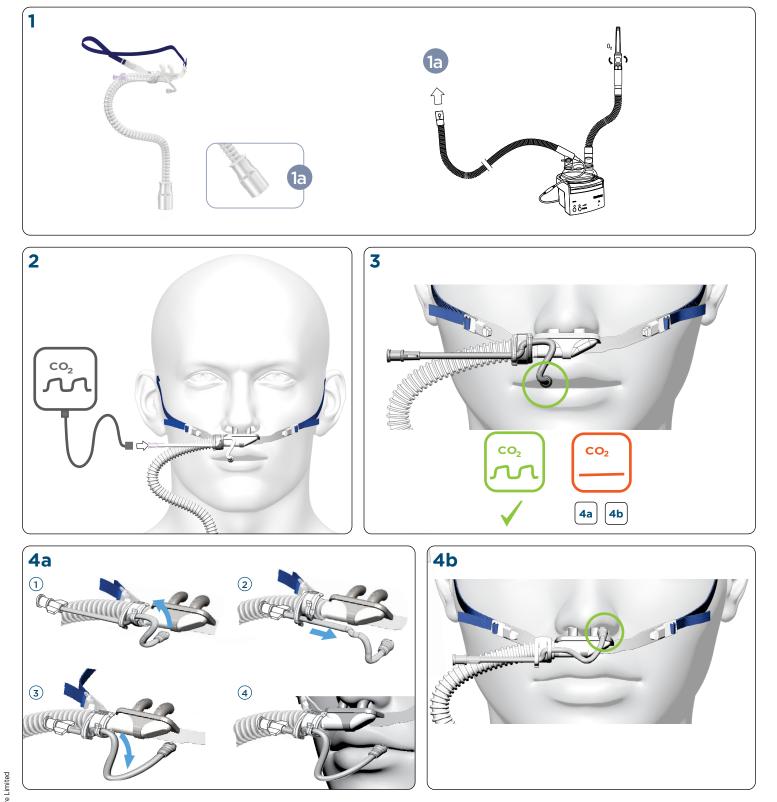

Rx only (2)

Fisher & Paykel Healthcare Limited. For patent information, see www.fphcare.com/ip

Optiflow

TRACE

F&

China 代理人/售后服务机构:费雪派克医疗保健(广州)有限公司,广州高新技术产业开发区科学域科丰路31号G12栋301号 电话: +86 20 32053486 传真: +86 20 32052132 Denmark Tel: +45 70 26 37 70 Fax: +46 83 66 310 Finland Tel: +358 9 251 66 123 Fax: +46 83 66 310 France Tel: +33 1 6446 5201 Finland Tel: +358 9 251 66 123 Fax: +46 83 66 310 France Tel: +33 1 6446 5201 Fax: +33 1 6446 5221 Email: c.s@fphcare.fr Germany Fisher & Paykel Healthcare GmbH, Deutschland, Österreich, Schweiz, Wiesenstrasse 49, D 73614 Schorndorf, Germany Tel: +49 7181 98599 0 Fax: +49 7181 98599 66 Hong Kong Tel: +852 2116 0032 Fax: +852 2116 0085 India Tel: +91 80 2309 6400 Ireland Tel: 1800 409 011 Fax: +44 1628 626 146 Italy Tel: +39 06 7839 2939 Fax: +39 06 7814 7709 Japan Tel: +81 3 5117 7110 Fax: +81 3 5117 7115 Korea Tel: +82 2 6205 6900 Fax: +82 2 6309 6901 Mexico Tel: +52 55 9130 1626 Norway Tel: +47 21 60 13 53 Fax: +47 22 99 60 10 Poland Tel: 00800 49 119 77 Fax: 00800 49 119 78 Russia Tel. and Fax: +7 495 782 2150 Spain Tel: +34 902 013 346 Fax: +34 902 013 379 Sweden Tel: +46 8 564 76 680 Fax: +46 8 36 6310 Switzerland Tel: 0800 83 47 63 Fax: 0800 83 47 54 Taiwan Tel: +886 2 87511739 Fax: +886 2 8751 5625 Turkey İthalatçı Firma: Fisher Paykel Sağlık Ürünleri Ticaret Limited Şirketi, İletişim Bilgileri: Ostim Mahallesi 1249. Cadde No:6, Yenimahalle, Ankara, Türkiye 06374, Tel: +90 312 354 34 12 Fax: +90 312 354 31 01 UK Fisher & Paykel Healthcare Ltd, Unit 16, Cordwallis Park, Clivemont Road, Maidenhead, Berkshire SL6 7BU, UK Tel: 0800 132 189 Fax: +44 1628 626 146 USA/Canada Tel: 1800 446 3908 or +1 949 453 4000 Fax: +1 949 453 4001 Fax: +1 949 Fax: +1 94 Fisher & Paykel Healthcare, F&P, Optiflow and Optiflow TRACE are trademarks of

Manufacturer and Fisher & Paykel Healthcare Ltd, 15 Maurice Paykel Place, East Tamaki, Auckland 2013, PO Box 14 348 Panmure, Auckland 1741, New Zealand Tel: +64 9 574 0100 Fax: +64 9 574 0158 Email: info@fphcare.co.nz Web: www.fphcare.com Importer/ Distributor Australia (Sponsor) Fisher & Paykel Healthcare Pty Ltd, 19-31 King Street, Nunawading, Melbourne, Victoria 3131. Tel: +61 3 9871 4900 Fax: +61 3 9871 4998 Austria Tel: 0800 29 31 23 Fax: 0800 29 31 22 Benelux Tel: +31 40 216 3555 Fax: +31 40 216 3555 Brazil Fisher & Paykel do Brasil, Rua Sampaio Viana, 277 cj 21, Paraíso, 04004-000, São Paulo – SP, Brazil Tel: +55 11 2548 7002



English (en)

F&P Optiflow^M Nasal Interface with CO₂ Sampling AA030J S/M/L (Small / Medium / Large)

Indications for Use

This product is a single-patient-use device that delivers respiratory gases to adult patients in hospitals.

This product is indicated for the delivery of Nasal High Flow (NHF) and Low Flow Oxygen to spontaneously breathing patients by appropriately qualified healthcare professionals.

Qualitative carbon dioxide sampling can be used at nasal cannula flow rates from 5 to 50 L/min in operating and procedure rooms

This product can be used for pre-oxygenation and short-term supplemental oxygenation (up to 10 minutes) during intubation in operating rooms under the direction of a physician anesthesiologist. This product is not indicated for apneic ventilation

∠L Contraindications

The following are contraindicated for this product. Failure to comply can lead to serious injury or death

- Do not use this product with electrosurgery or electrocautery devices on the head or neck. Refer to the FIRE DANGER information.
- Do not use where Continuous Positive Airway Pressure (CPAP) is contraindicated (e.g. pneumothorax, bullous lung disease, craniofacial trauma or surgery, airway foreign body, unstable hemodynamics). Do not use where nasal interfaces are contraindicated (e.g. nasal
- obstruction, nasal trauma). Do not administer inhalants, anesthetic gases or vapors through this
- product.
- Do not use an anesthesia mask while this product is in use. Attempting to seal the mask over the product may result in excessive delivered pressure.

I FIRE DANGER

This product is an open oxygen delivery system. Open oxygen delivery can increase the risk of a surgical fire occurring, causing serious injury or death. Extreme care must be taken. The following is advised.

- Contraindication: Do not use this product with electrosurgery or electrocautery devices on the head or neck.
- Warning: Do not use this product where ignition sources and fuel are present. Use with ignition sources and fuel present completes the fire triangle, increasing the risk of fire. The following steps should be taken to reduce the risk of a surgical fire occurring:
- Evaluate the oxygen needs of each individual patient and use the minimum supplementation required. Oxygen titration with a gas mixer should be considered.
- Ensure the use of this product is communicated to all operating room personnel when conducting the procedure's fire risk assessment.
- Ensure gas flow from the product is not in or near the surgical field AND the surgical site is free of all potential fuel sources, including alcohol skin preparation, gauze, sponges and drapes BEFORE potential ignition sources, such as electrosurgery, electrocautery or laser devices are used.
- Ensure the room is adequately ventilated as oxygen may accumulate over time.
- Ensure gas flow from the product does not pool under drapes. Follow the instructions for use of all surgical devices, including electrosurgery, electrocautery and laser devices, regarding oxygen delivery.

Specifications

- Product is for single-patient-use only for a maximum period of 24 hours. Operating flow range: 5 – 70 L/min.
- **Caution:** Flow above 15 L/min must be humidified. Failure to humidify can cause nasal and mucosal dryness and discomfort.
- Carbon dioxide sampling only functions within the operating flow range of 5 50 L/min. Note: Higher flow rates may dilute the qualitative trace below the minimum level of detection of the carbon dioxide analyzer.
- Use only with compatible Fisher & Paykel Healthcare (F&P) humidifiers and accessories. **Warning:** Incompatible humidifiers and accessories which are used with this product may impair the performance of this product or compromise safety (including potentially causing serious injury or death). Ambient operating temperature: 18 – 26 °C.
- This product is not intended to be sterile.
- This product was not made with natural rubber latex.

Setup Instructions

Connect product to compatible system.

- 2. Connect carbon dioxide sampling tube to carbon dioxide sampling line and analyzer. **Warning:** Luer connection for gas sampling only. Do not connect to intravenous fluids or any other device.
- Ensure the recommended water prevention system for the carbon dioxide analyzer is installed. **Note:** Failure to do so may cause condensed water, secretions and/or contaminates to impair the analyzer's performance.
- Fit product to patient.
- Ensure product is sized correctly. Warning: Do not create a seal in the nares as this may present a risk of barotrauma or gastric insufflation.
- Ensure all connections are secure.
- Ensure a carbon dioxide trace is present on the carbon dioxide analyzer to confirm carbon dioxide sampling tube function and location. If a trace cannot be obtained refer to the troubleshooting instructions.
- Set flow rate and ensure gas flow is exiting the prongs

Troubleshooting Instructions

- If carbon dioxide sampling tube is not long enough to reach your desired sampling location:
 - Unwind tube from plastic clip.
 - (2) Pull tube to desired length.
 - (3) Rewind tube onto plastic clip. (4) Bend tube to desired sampling location.
- 4b. If a qualitative sample cannot be obtained, bend carbon dioxide sampling tube towards patient's nose.

Operating Instructions

During use of this product, patients may desaturate rapidly with little warning and / or their arterial carbon dioxide may increase. Warning: Failure to follow the practices described below and appropriately address these situations can lead to serious injury or death.

- Use appropriate patient monitoring, including oxygen saturation and carbon dioxide monitoring.
- For spontaneously breathing patients, this product provides a qualitative carbon dioxide trace only. Alternative monitoring methods should be used where indicated.
- Ensure availability of appropriate measures to respond to a rapid desaturation and/or carbon dioxide increase
- Remove the product from the patient should an anesthesia mask based therapy be required.
- In addition, during apnea, the following practices must also be adhered to:
- Maintain a patent airway at all times (e.g. by applying jaw thrust).
- Do not use product where a secure airway is indicated. This may increase the risk of pulmonary aspiration

/!\ Warnings

- Do not occlude the patient's mouth or nose with hands or any devices. Occlusion can lead to excessive delivered pressure and cause serious injury or death.
- This product delivers oxygen. Oxygen titration with a gas mixer should be considered where hyperoxia is a concern.

∠L Cautions

- Do not use product where quantitative exhaled gas measurements are indicated and a second and the second and the second and the second are indicated as a second and the second as a second a arterial blood gas test should be considered.
- Do not reuse product on multiple patients. Reuse may result in transmission of infectious substances, interruption to treatment, serious injury or death.
- Do not use product if damaged. Product damage may result in impaired performance.
- Do not remove the clip holding the carbon dioxide sampling tube from the product as this may damage the product or cause leaks in the system and reduce the gas supply to the patient.
- Do not soak, wash or sterilize the product. Avoid contact with chemicals, cleaning agents and hand sanitizers as this may damage the product or cause leaks in the system and reduce the gas supply/ humidity to the patient.
- Do not stretch or crush tubing, as this may occlude the gas flow or damage the tubing, resulting in impaired performance or gas leaks Do not use product if packaging is not sealed or it is contaminated, as this presents a risk of infection.
- Do not use product in or near a magnetic resonance imaging (MRI) scanner
- Dispose of product safely in accordance with standard hospital procedure. User may be exposed to respiratory tract fluids during disposal.

Notes

- Refer to the instructions for use of compatible Fisher & Pavkel Healthcare humidifiers and accessories for additional warnings, cautions, contraindications, and system information.
- If a serious incident has occurred while using this product, please contact your local Fisher & Paykel Healthcare representative and Competent Authority.