




The F&P Eson's three key components, the **RollFit™ Seal**, **ErgoFit Headgear** and **Easy Frame** work in harmony. These three key components have been designed to deliver the comfort, seal and easy use that Fisher & Paykel Healthcare masks are known for.

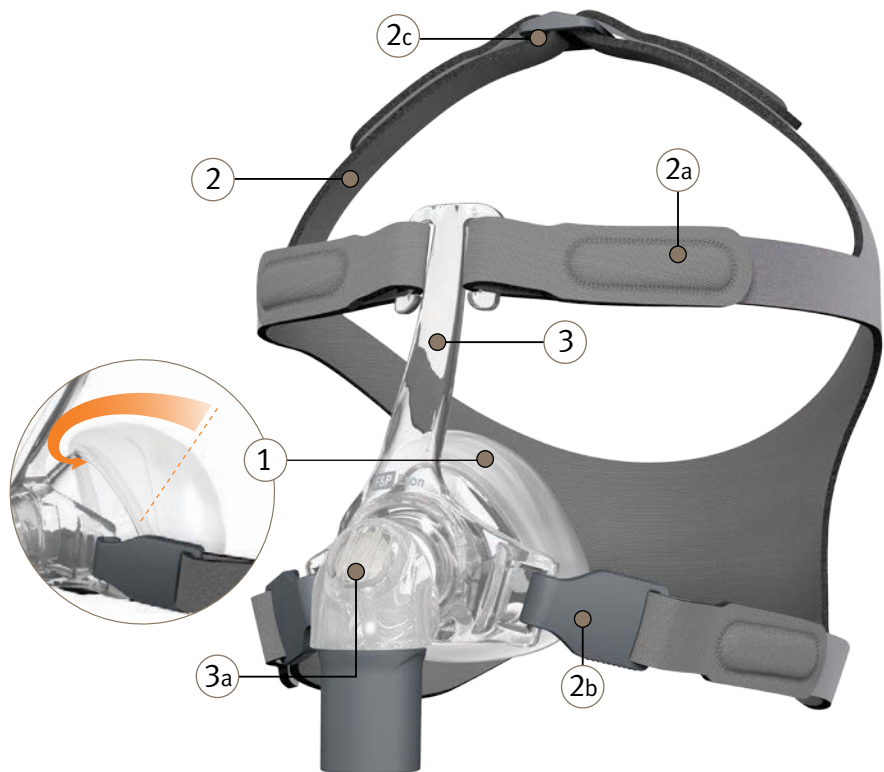

### 1. RollFit Seal

As the name suggests, the one-piece seal 'rolls' back and forth on the bridge of the nose to adjust automatically. In doing so, the **RollFit Technology** reduces pressure on the bridge of the nose.

The **RollFit Seal** enhances sealing performance and stability.

### 2. ErgoFit Headgear

This breathable, ergonomically designed headgear has been made to self-locate high on the rear of the head, allowing for head movement (sideways, up and down).

The **ErgoFit Headgear** also incorporates embedded **VELCRO®\*** **Tabs** (2a) with **Easy-Clip hooks** (2b) for fast fitting and removal.

The adjustable **Crown Strap** (2c) enhances mask stability.

\* VELCRO® is a registered trademark of Velcro Industries B.V.

### 3. Easy Frame

This low-profile frame is stable, durable and small, and ensures a clear line of sight.

The one frame fits all three seal sizes and has an **Easy-Clip** frame attachment to aid assembly after cleaning.

The **Easy Frame** comes with the **Q Cover** (3a) for quiet use<sup>1</sup> and minimal partner draft.

The **Ball and Socket Elbow** rotates for freedom of movement while reducing drag from the CPAP tube.

Mask Reorder Codes		Spare Parts Reorder Codes			
<b>400449</b>	F&P Eson Mask – Small	<b>400HC567</b>	F&P Eson Headgear – Small	<b>400HC228</b>	F&P Eson Diffuser (10pk) and Cover
<b>400450</b>	F&P Eson Mask – Medium	<b>400HC568</b>	F&P Eson Headgear – Medium/Large	<b>400HC570</b>	F&P Eson Swivel
<b>400451</b>	F&P Eson Mask – Large	<b>400HC569</b>	F&P Eson Headgear – Clips and Buckle	<b>400HC571</b>	F&P Eson Mask No H/G – Small
		<b>400HC126</b>	F&P Eson Seal – Small	<b>400HC572</b>	F&P Eson Mask No H/G – Medium
		<b>400HC127</b>	F&P Eson Seal – Medium	<b>400HC573</b>	F&P Eson Mask No H/G – Large
		<b>400HC128</b>	F&P Eson Seal – Large		

1. The Sound Pressure Level of the mask with the Q Cover is 22 dBA, with uncertainty 2.5 dBA.