

Patient summary report - 16 Dec 2016 - 14 Jan 2017

PATIENT DETAILS

Name:	Report Demo	Phone:	0123456789
Patient ID:	1234	Email:	sample@example.com
DOB:	07 Jan 1960	Gender:	Other

Physician:	Example GP	Phone:	0123456789
------------	------------	--------	------------

Avg. AHI:	Avg. Usage (h):	Compliance*:	* Last 30 Days
5.8	5.6	87 %	% of nights CPAP used ≥ 4 hours from 16 Dec 2016

PATIENT DEVICE

AS AT 14/01/2017

Device:	SleepStyle Auto	Humidity:	5.0
Serial No:	1111111111	Expiratory Relief:	2 cmH ₂ O
Mode:	Auto	SensAwake:	On
Min Set Pressure:	6.0 cmH ₂ O	SensAwake Pressure:	6.0 cmH ₂ O
Max Set Pressure:	13.0 cmH ₂ O		



COMPLIANCE

30 DAY PERIOD

16 Dec 2016 - 14 Jan 2017

FROM SETUP DATE

19 Jul 2016 - 14 Jan 2017

Total Days:

Used: 27 (90 %) 168 (93 %)

Used ≥ 4 hrs: 26 (87 %) 160 (89 %)

Not Used: 3 (10 %) 12 (7 %)

Average Usage (Days on Therapy):

6.3 h 7.0 h

Average Usage (All Days):

5.6 h 6.5 h

Total Run Time:

169.0 h 1188.2 h

EFFICACY**

AHI: 5.8

Apnea Length: 13 s

CAI: 0.5

Pressure: 7 cmH₂O90th % Pressure: 10 cmH₂O

At Max Pressure: 4 %

System Leak: 42 l/min

90th % Leak: 62 l/min

Time With Excessive Leak: 13 %

No. of SensAwake Detections: 17 per night

Time with SensAwake on: 100 %

ENVIRONMENT

Ambient Temperature: 24 °C

USAGE

hours


Average : 5.6

Days Used : 90 %

Total Run Time

Total Usage

Usage < 4 h

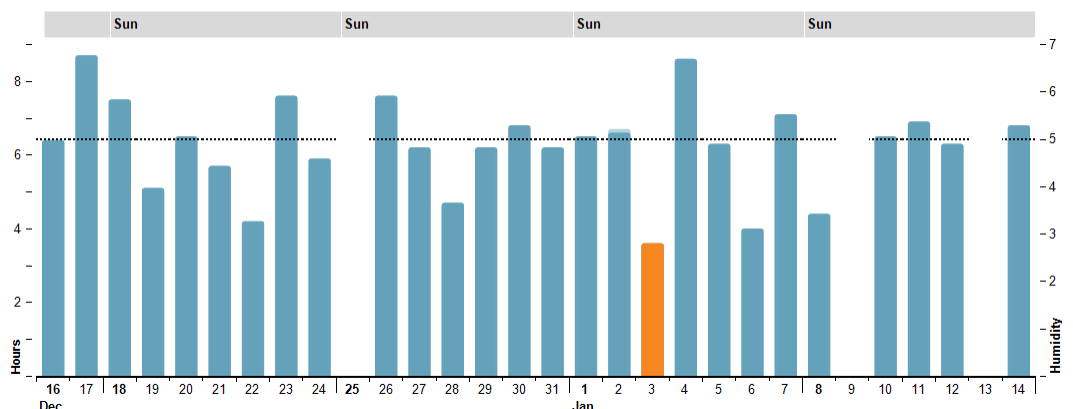

Patient summary report - 16 Dec 2016 - 14 Jan 2017

PATIENT DETAILS


Name: Report Demo Phone: 0123456789
Patient ID: 1234 Email: sample@example.com
DOB: 07 Jan 1960 Gender: Other

AHI

Events per hour

AHI : 5.8
CAI : 0.5
OAI : 1.1
HI : 4.2

Total AHI
AHI

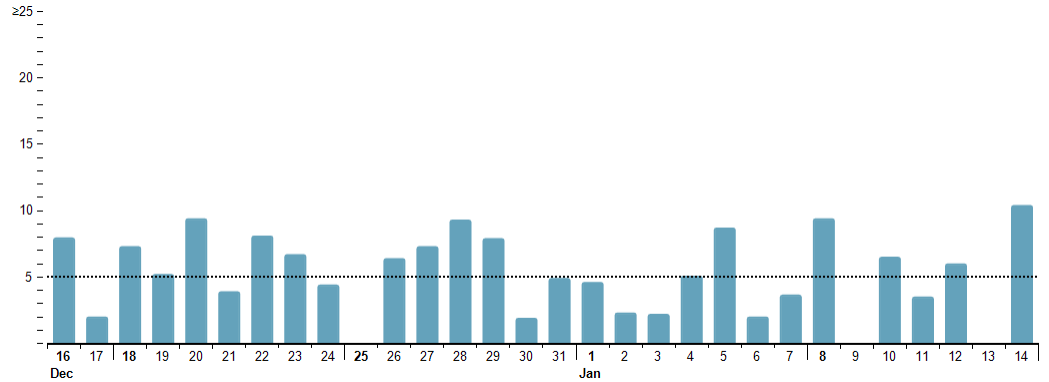


PRESSURE

cmH₂O

Average : 7.0
90th Percentile : 10.0

Average
90th Percentile
Max Pressure Reached
Min / Max

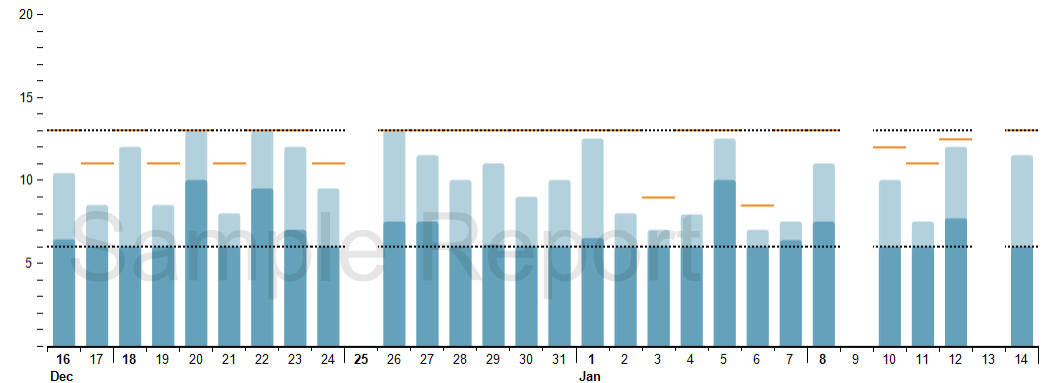


SYSTEM LEAK

l/min

Average : 42
90th Percentile : 62

Average
90th Percentile
High Leak

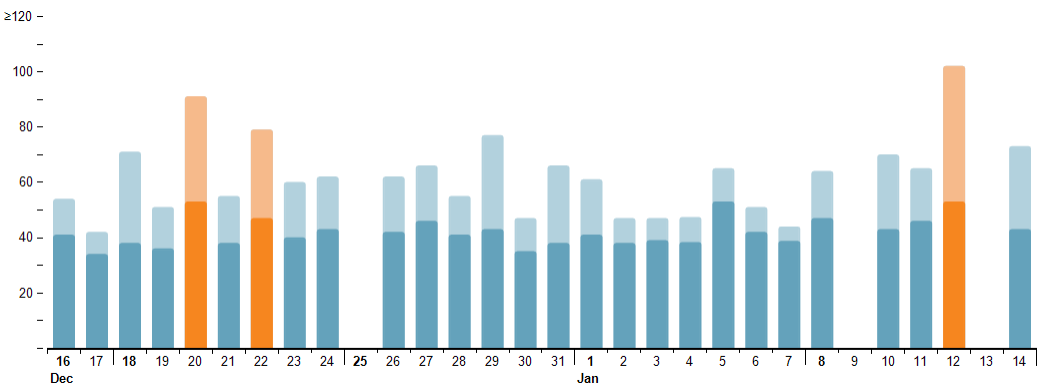


USAGE

hours

Average : 5.6
Days Used : 90 %

Total Run Time
Total Usage
Usage < 4 h

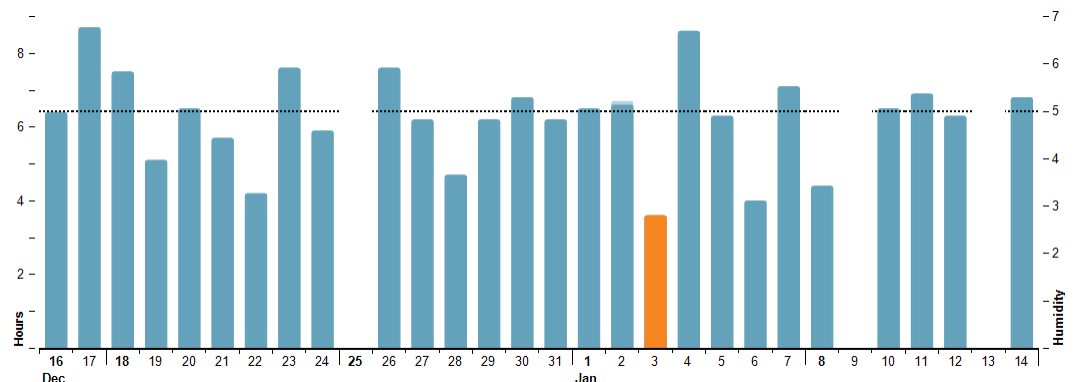


Patient summary report - 16 Dec 2016 - 14 Jan 2017

COMPLIANCE SUMMARY

Name:	Report Demo	Phone:	0123456789
Patient ID:	1234	Email:	sample@example.com
DOB:	07 Jan 1960	Gender:	Other

Prescription history

PMin: 6.0 PMax: 13.0
 H: 5.0 SA: 6.0 ER: 2.0

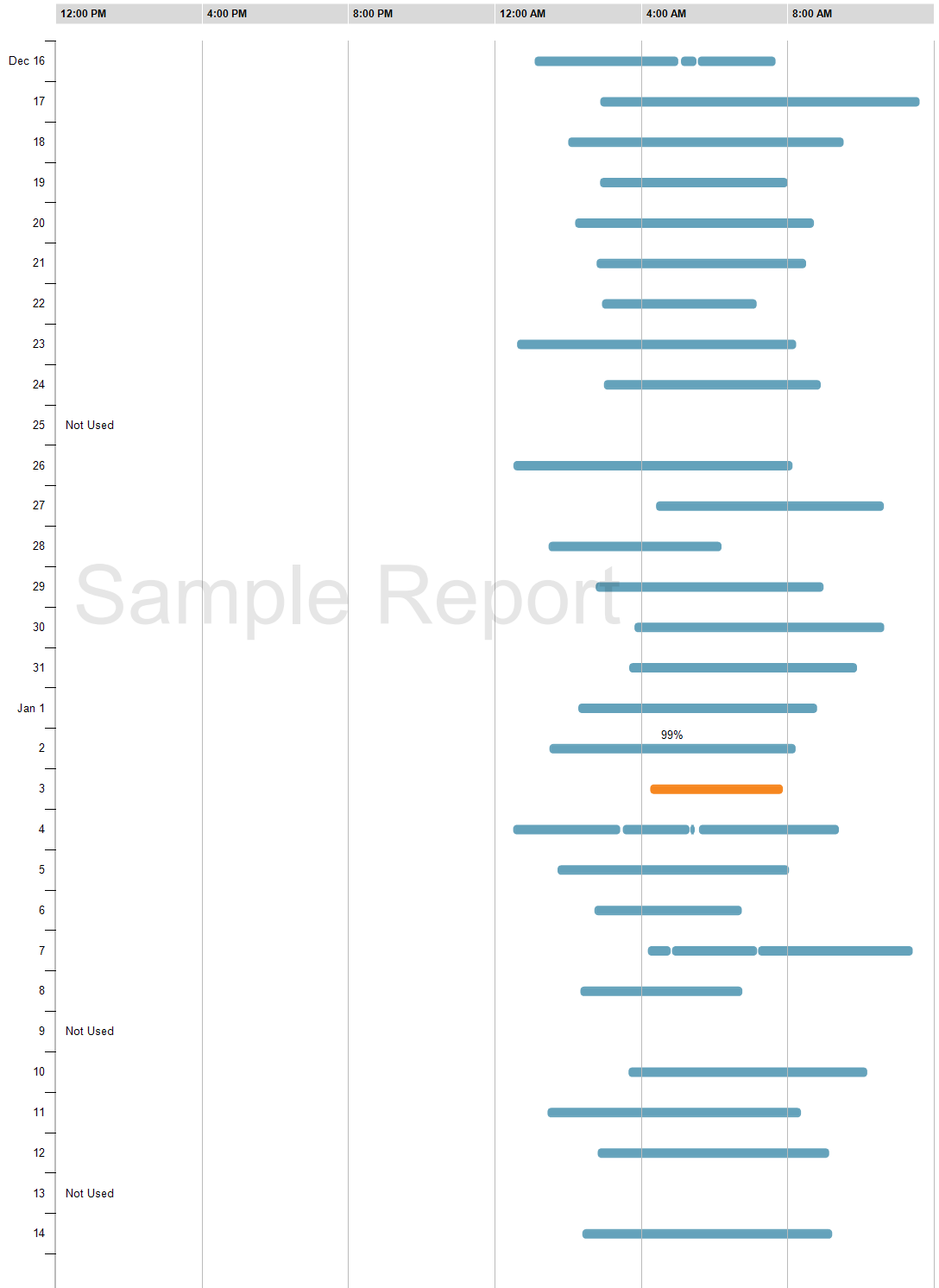

■ Compliance < 4 hours
 ■ Compliance < 4 hours and < 100% usage
 ■ 100% usage
 ■ < 100% usage
 82%
 82%

PMin: Min Pressure PMax: Max Pressure H: Humidity SA: SensAwake ER: Expiratory Relief
 Prescription Change

Patient summary report - 16 Dec 2016 - 14 Jan 2017

08 Jan 2017

Sunday 2:20 AM - 6:44 AM

Name: Report Demo

Patient ID: 1234

DOB: 07 Jan 1960

Phone: 0123456789


Email: sample@example.com

Gender: Other

RUN TIME

Total Run Time : 4.4 h

Percentage Usage : 100 %

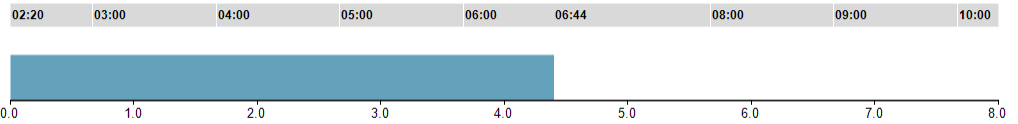

EVENTS

AHI : 9.4

Total Obstructive Apnea : 11

OAI : 2.5

Total Hypopneas : 28

HI : 6.4

Total Central Apneas : 2


CAI : 0.5

Total Flow Limitation : 18

Flow Limitation /Hr : 4.1

Total SensAwake : 12

SensAwake /Hr : 3.0

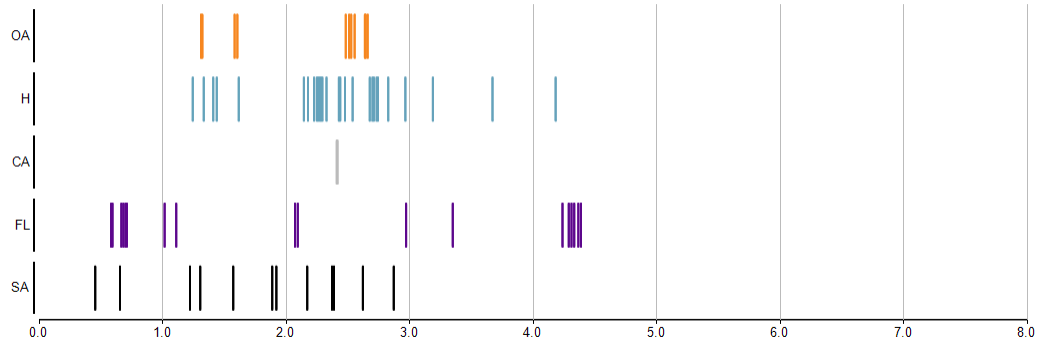

PRESSURE


cmH₂O

Median : 7.5

90th Percentile : 11.0

/// Ramp Period

..... Min / Max

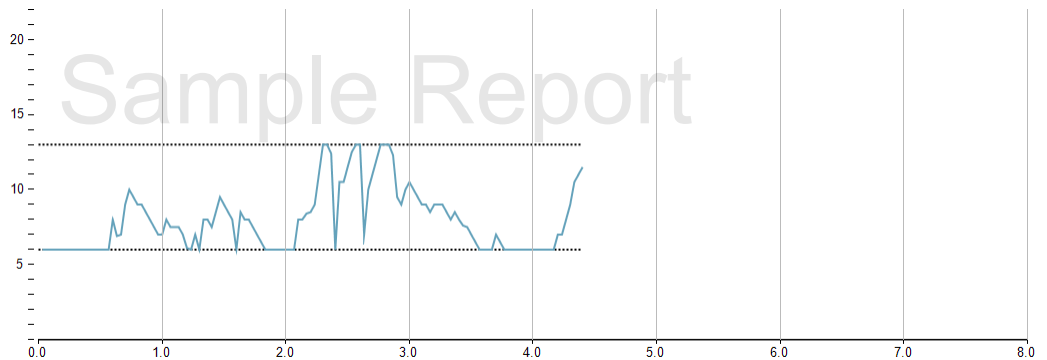


SYSTEM LEAK

l/min

Median : 47

90th Percentile : 64

Time with Excessive Leak : 17 %

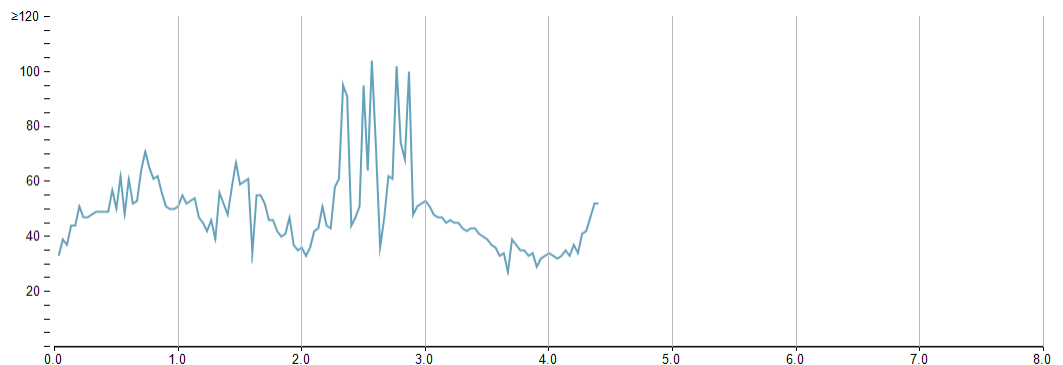

TIME AND EVENTS AT PRESSURE

P (cmH ₂ O)	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20
T (min)	0	0	94	30	44	40	16	14	10	16	0	0	0	0	0	0	0
OA	0	0	1	0	2	1	3	0	1	3	0	0	0	0	0	0	0
H	0	0	0	3	2	5	2	6	3	7	0	0	0	0	0	0	0
CA	0	0	0	0	0	0	0	2	0	0	0	0	0	0	0	0	0
FL	0	0	0	1	7	4	1	4	1	0	0	0	0	0	0	0	0


Patient summary report - 16 Dec 2016 - 14 Jan 2017

10 Jan 2017

Tuesday 3:39 AM - 10:09 AM


Name: Report Demo
Patient ID: 1234
DOB: 07 Jan 1960Phone: 0123456789
Email: sample@example.com
Gender: Other

RUN TIME

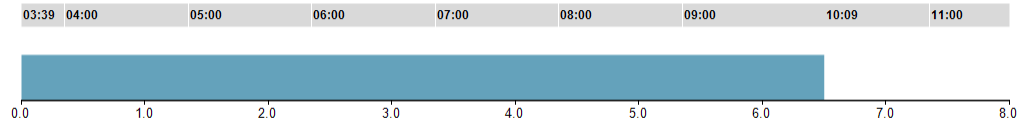Total Run Time : 6.5 h
Percentage Usage : 100 %

EVENTS

AHI : 6.5

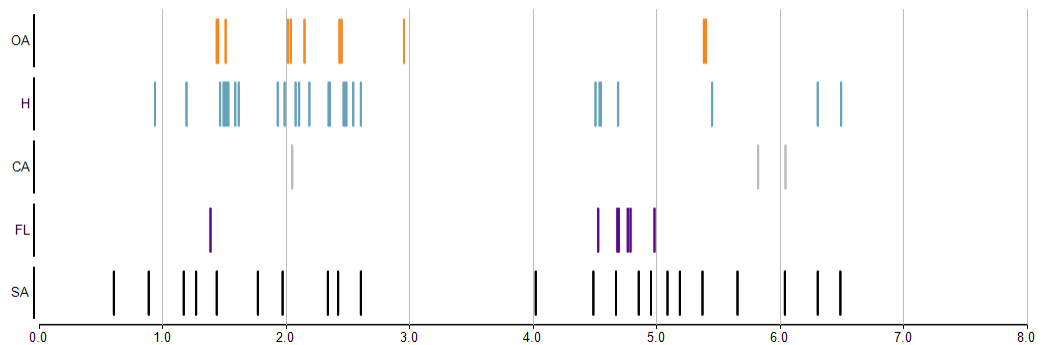Total Obstructive Apnea : 11
OAI : 1.7Total Hypopneas : 28
HI : 4.3Total Central Apneas : 3
CAI : 0.5Total Flow Limitation : 7
Flow Limitation /Hr : 1.1Total SensAwake : 22
SensAwake /Hr : 3.4

PRESSURE


cmH₂O

Median : 6.0

90th Percentile : 10.0

Ramp Period

Min / Max

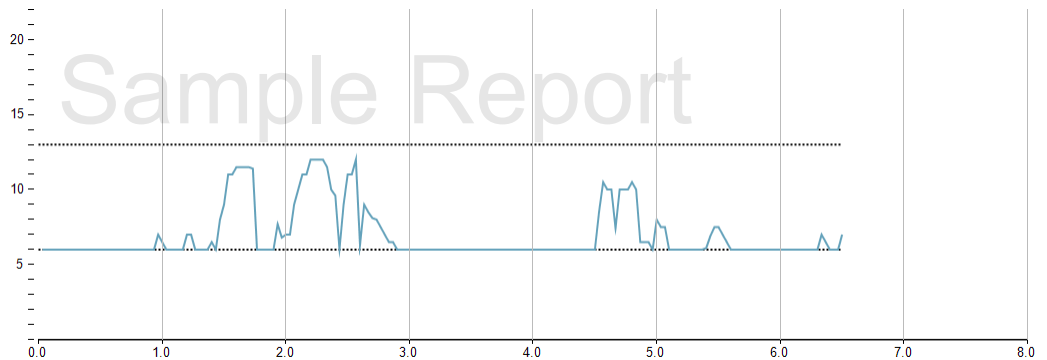


SYSTEM LEAK

l/min

Median : 43

90th Percentile : 70

Time with Excessive Leak : 14 %

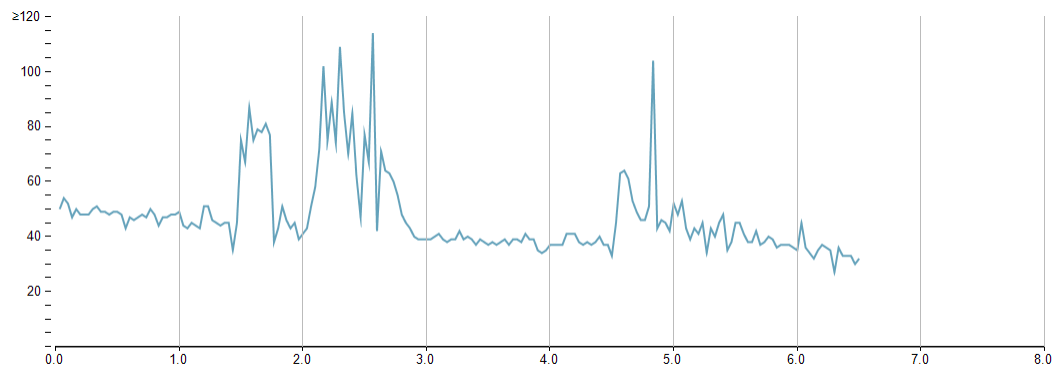

TIME AND EVENTS AT PRESSURE

P (cmH ₂ O)	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20
T (min)	0	0	260	42	20	12	18	18	20	0	0	0	0	0	0	0	0
OA	0	0	4	1	2	2	0	2	0	0	0	0	0	0	0	0	0
H	0	0	0	5	3	4	4	8	4	0	0	0	0	0	0	0	0
CA	0	0	2	0	0	1	0	0	0	0	0	0	0	0	0	0	0
FL	0	0	0	1	1	1	3	1	0	0	0	0	0	0	0	0	0

Patient summary report - 16 Dec 2016 - 14 Jan 2017

11 Jan 2017

Wednesday 1:26 AM - 8:20 AM

Name:

Report Demo

Phone:

0123456789

Patient ID:

1234

Email:

sample@example.com

DOB:

07 Jan 1960


Gender:

Other

RUN TIME

Total Run Time : 6.9 h

Percentage Usage : 100 %

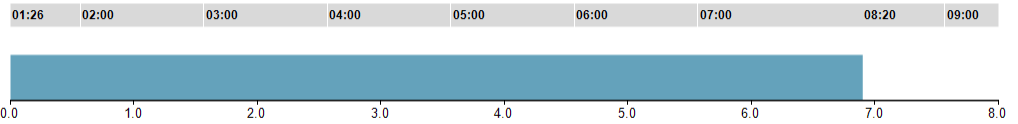

EVENTS

AHI : 3.5

Total Obstructive Apnea : 6

OAI : 0.9

Total Hypopneas : 13

HI : 1.9

Total Central Apneas : 5


CAI : 0.7

Total Flow Limitation : 7

Flow Limitation /Hr : 1.0

Total SensAwake : 15

SensAwake /Hr : 2.2

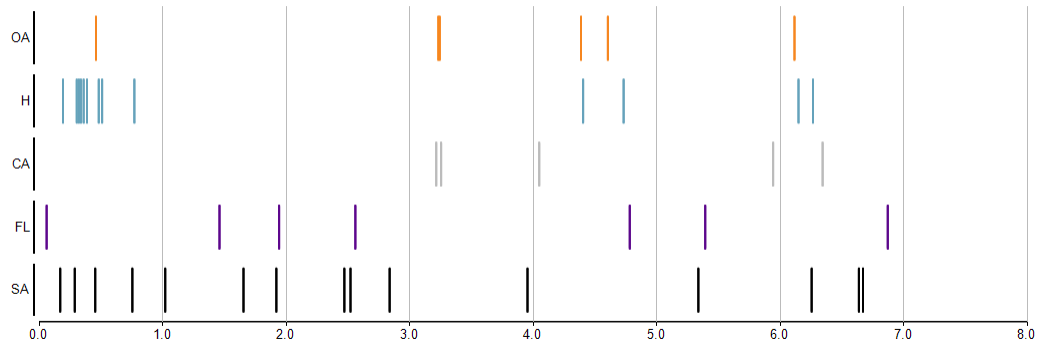

PRESSURE


cmH₂O

Median : 6.0

90th Percentile : 7.5

Ramp Period

Min / Max

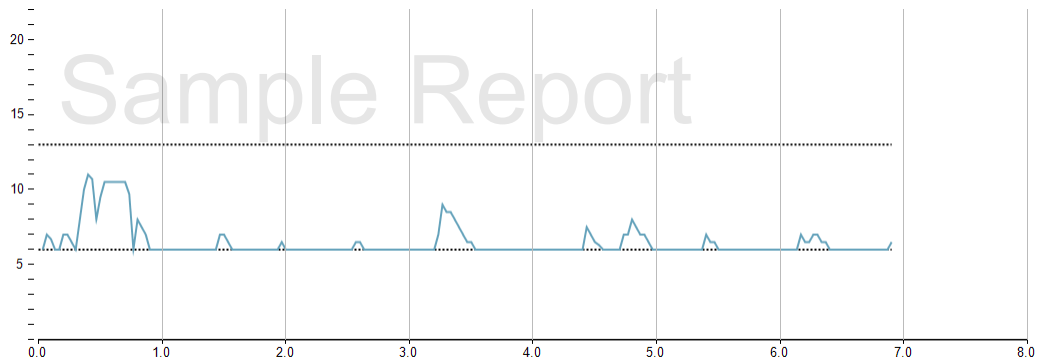


SYSTEM LEAK

l/min

Median : 46

90th Percentile : 65

Time with Excessive Leak : 13 %

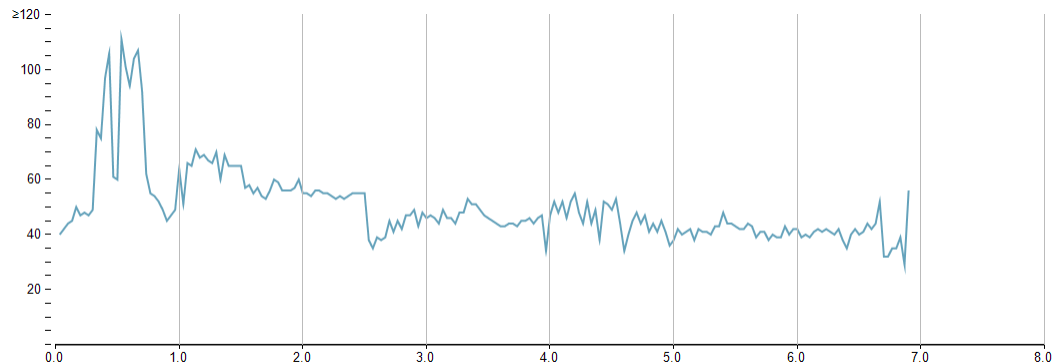

TIME AND EVENTS AT PRESSURE

P (cmH ₂ O)	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20
T (min)	0	0	300	68	18	6	6	16	0	0	0	0	0	0	0	0	0
OA	0	0	3	1	1	1	0	0	0	0	0	0	0	0	0	0	0
H	0	0	0	4	4	0	3	2	0	0	0	0	0	0	0	0	0
CA	0	0	2	2	0	1	0	0	0	0	0	0	0	0	0	0	0
FL	0	0	0	6	1	0	0	0	0	0	0	0	0	0	0	0	0

Patient summary report - 16 Dec 2016 - 14 Jan 2017

12 Jan 2017


Thursday 2:48 AM - 9:06 AM

Name: Report Demo
Patient ID: 1234
DOB: 07 Jan 1960

Phone: 0123456789
Email: sample@example.com
Gender: Other


RUN TIME

Total Run Time : 6.3 h
Percentage Usage : 100 %

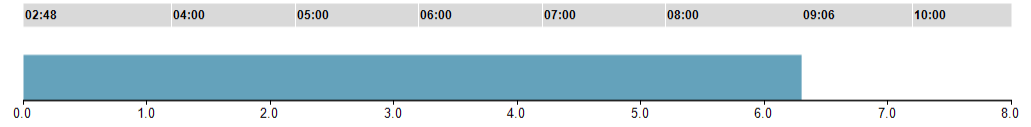

EVENTS


AHI : 6.0
Total Obstructive Apnea : 9
OAI : 1.4
Total Hypopneas : 26
HI : 4.1
Total Central Apneas : 3
CAI : 0.5
Total Flow Limitation : 9
Flow Limitation /Hr : 1.4
Total SensAwake : 8
SensAwake /Hr : 1.3

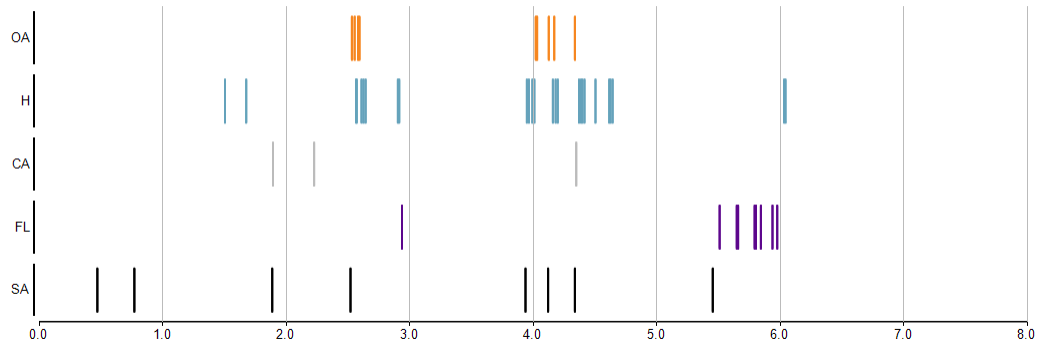

PRESSURE

cmH₂O


Median : 7.8
90th Percentile : 12.0
Ramp Period
Min / Max

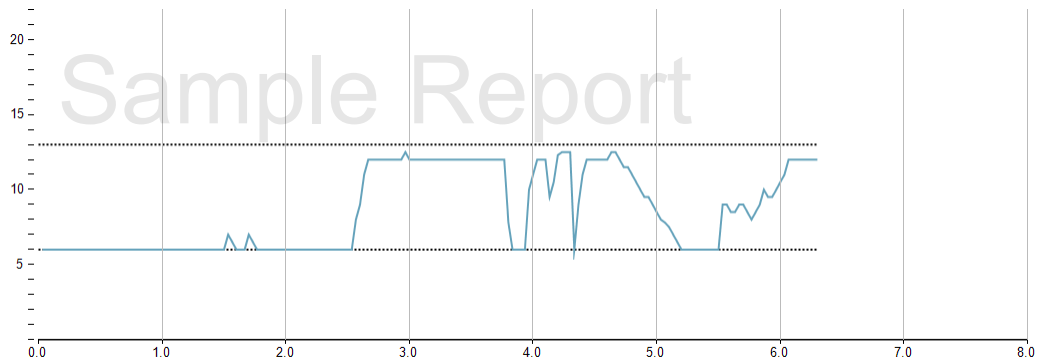

SYSTEM LEAK

l/min

Median : 53
90th Percentile : 102
Time with Excessive Leak : 38 %

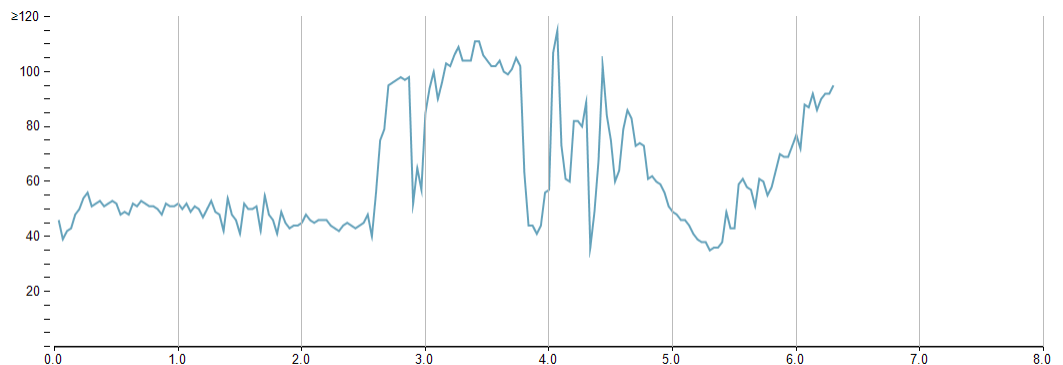

TIME AND EVENTS AT PRESSURE

P (cmH ₂ O)	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20
T (min)	0	0	174	12	12	26	18	16	108	12	0	0	0	0	0	0	0
OA	0	0	1	0	1	3	1	0	3	0	0	0	0	0	0	0	0
H	0	0	0	2	1	1	2	8	9	3	0	0	0	0	0	0	0
CA	0	0	2	0	0	1	0	0	0	0	0	0	0	0	0	0	0
FL	0	0	0	0	0	5	2	1	0	1	0	0	0	0	0	0	0

Patient summary report - 16 Dec 2016 - 14 Jan 2017

14 Jan 2017

Saturday 2:23 AM - 9:11 AM

Name: Report Demo

Patient ID: 1234

DOB: 07 Jan 1960

Phone: 0123456789


Email: sample@example.com

Gender: Other

RUN TIME

Total Run Time : 6.8 h

Percentage Usage : 100 %

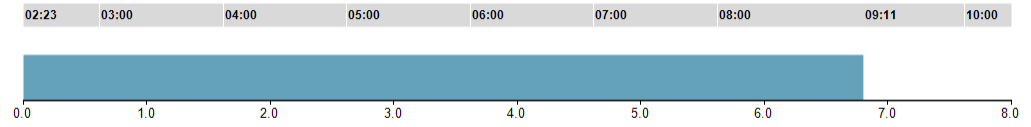

EVENTS

AHI : 10.4

Total Obstructive Apnea : 31

OAI : 4.6

Total Hypopneas : 24

HI : 3.6

Total Central Apneas : 15


CAI : 2.2

Total Flow Limitation : 0

Flow Limitation /Hr : 0.0

Total SensAwake : 21

SensAwake /Hr : 3.1

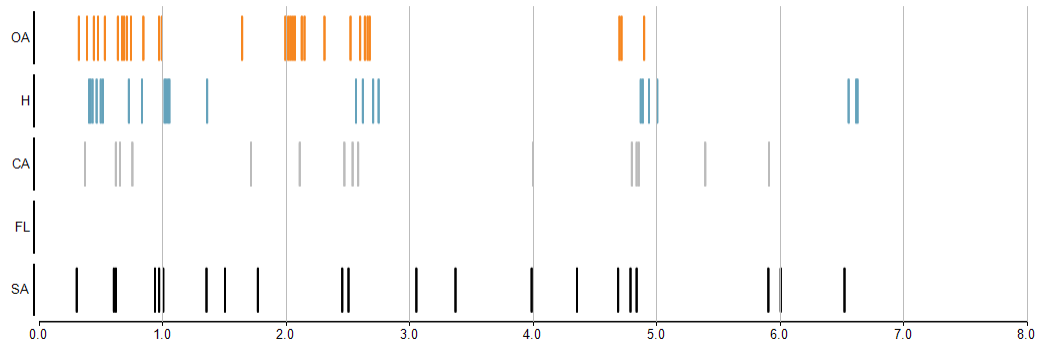

PRESSURE


cmH₂O

Median : 6.0

90th Percentile : 11.5

Ramp Period

Min / Max

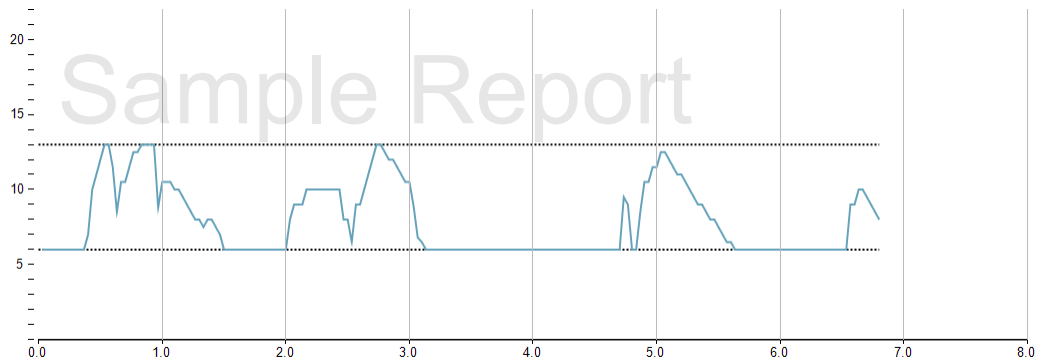


SYSTEM LEAK

l/min

Median : 43

90th Percentile : 73

Time with Excessive Leak : 23 %

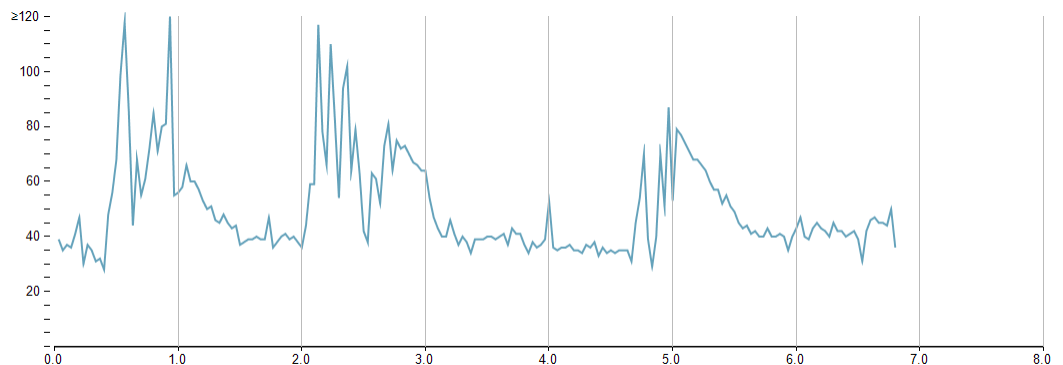

TIME AND EVENTS AT PRESSURE

P (cmH ₂ O)	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20
T (min)	0	0	210	16	26	38	40	30	22	26	0	0	0	0	0	0	0
OA	0	0	4	2	2	5	3	9	3	3	0	0	0	0	0	0	0
H	0	0	0	0	1	2	6	7	3	5	0	0	0	0	0	0	0
CA	0	0	4	1	1	7	0	1	0	1	0	0	0	0	0	0	0
FL	0	0	0	0	0	0	0	0	0	0	0	0	0	0	0	0	0

27 DAY ANALYSIS

Name: Report Demo
Patient ID: 1234
DOB: 07 Jan 1960
Phone: 0123456789
Email: sample@example.com
Gender: Other

Date	Mode	Pressure (cmH ₂ O)	Humidity	Ramps	Total Run Time	Usage Time (h)	AHI	CAI	Number of SensAwake Detections	SensAwake Pressure (cmH ₂ O)	Expiratory Relief	90th Percentile Pressure (cmH ₂ O)	Max Pressure Reached (cmH ₂ O)	90th Percentile Leak (l/min)
16 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	6.4	6.4	8.0	0.2	22	6.0	2	10.5	13.0	54
17 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	8.7	8.7	2.0	0.3	22	6.0	2	8.5	11.0	42
18 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	7.5	7.5	7.3	0.4	25	6.0	2	12.0	13.0	71
19 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	5.1	5.1	5.2	0.0	16	6.0	2	8.5	11.0	51
20 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	6.5	6.5	9.4	0.6	14	6.0	2	13.0	13.0	91
21 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	5.7	5.7	3.9	0.2	18	6.0	2	8.0	11.0	55
22 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	4.2	4.2	8.1	1.2	12	6.0	2	13.0	13.0	79
23 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	7.6	7.6	6.7	0.3	18	6.0	2	12.0	13.0	60
24 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	5.9	5.9	4.4	0.3	14	6.0	2	9.5	11.0	62
26 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	7.6	7.6	6.4	0.1	21	6.0	2	13.0	13.0	62
27 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	6.2	6.2	7.3	0.0	16	6.0	2	11.5	13.0	66
28 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	4.7	4.7	9.3	0.4	17	6.0	2	10.0	13.0	55
29 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	6.2	6.2	7.9	0.6	28	6.0	2	11.0	13.0	77
30 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	6.8	6.8	1.9	0.1	12	6.0	2	9.0	13.0	47
31 Dec 2016	Auto	Min 6.0 Max 13.0	5.0	0	6.2	6.2	4.9	0.5	14	6.0	2	10.0	13.0	66
01 Jan 2017	Auto	Min 6.0 Max 13.0	5.0	0	6.5	6.5	4.6	0.6	15	6.0	2	12.5	13.0	61
02 Jan 2017	Auto	Min 6.0 Max 13.0	5.0	0	6.7	6.6	2.3	0.2	15	6.0	2	8.0	13.0	47
03 Jan 2017	Auto	Min 6.0 Max 13.0	5.0	0	3.6	3.6	2.2	0.0	10	6.0	2	7.0	9.0	47
04 Jan 2017	Auto	Min 6.0 Max 13.0	5.0	0	8.6	8.6	5.1	0.6	18	6.0	2	8.0	13.0	47
05 Jan 2017	Auto	Min 6.0 Max 13.0	5.0	0	6.3	6.3	8.7	0.3	16	6.0	2	12.5	13.0	65
06 Jan 2017	Auto	Min 6.0 Max 13.0	5.0	0	4.0	4.0	2.0	0.5	10	6.0	2	7.0	8.5	51
07 Jan 2017	Auto	Min 6.0 Max 13.0	5.0	0	7.1	7.1	3.6	0.8	17	6.0	2	7.5	13.0	44
08 Jan 2017	Auto	Min 6.0 Max 13.0	5.0	0	4.4	4.4	9.4	0.5	13	6.0	2	11.0	13.0	64
10 Jan 2017	Auto	Min 6.0 Max 13.0	5.0	0	6.5	6.5	6.5	0.5	22	6.0	2	10.0	12.0	70
11 Jan 2017	Auto	Min 6.0 Max 13.0	5.0	0	6.9	6.9	3.5	0.7	15	6.0	2	7.5	11.0	65
12 Jan 2017	Auto	Min 6.0 Max 13.0	5.0	0	6.3	6.3	6.0	0.5	8	6.0	2	12.0	12.5	102
14 Jan 2017	Auto	Min 6.0 Max 13.0	5.0	0	6.8	6.8	10.4	2.2	21	6.0	2	11.5	13.0	73

Note: The therapy data presented in this brochure was artificially generated and may not be representative or can omit real-life therapy data.