


Improving care for tracheostomy patients



Tracheostomy patients and their care providers no longer need to settle for the noise and sub-optimal humidification provided by compressors and cold-water nebulizers.

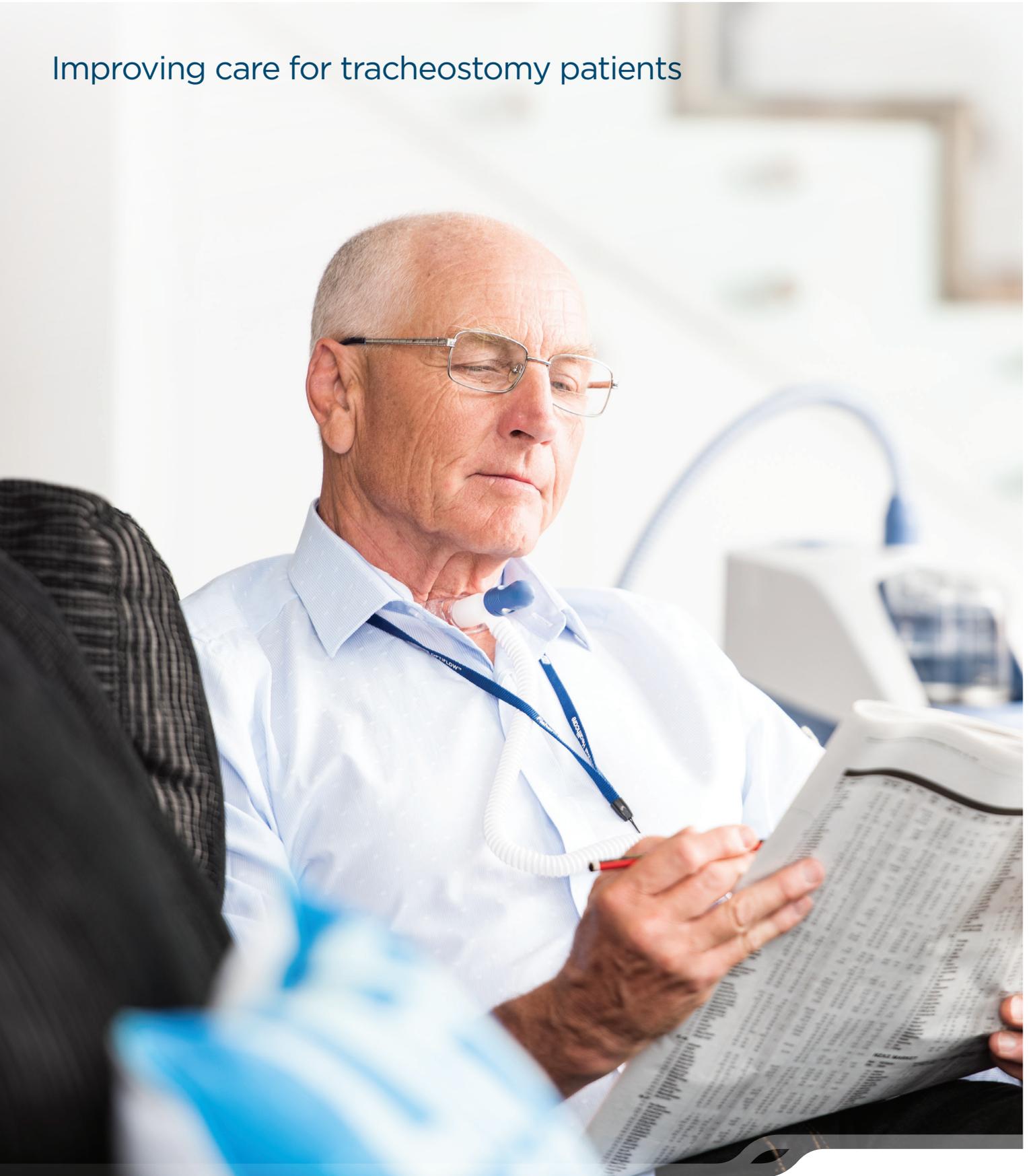

Optiflow™+ Tracheostomy interface

- Evaqua™ technology, which reduces the formation of mobile condensate in the patient interface
- Options include both direct-connect and tracheostomy collar interfaces

Supplementary oxygen (when required)

- FiO₂ 21% to 100%
- Inbuilt ultrasonic oxygen analyzer requires no calibration, service or replacement

Adjustable flow settings

- Integrated flow generator quietly delivers a wide flow range (10-60 L/min). No compressor or wall air supply required

Designed for simple setup, use and cleaning


- Helpful onscreen animations assist with setup and troubleshooting

Heated breathing tube (AirSpiral™)

- Dual spiral heater wires and unique integrated temperature sensor

Integrated temperature sensor


- No external probes, cables or adaptors required

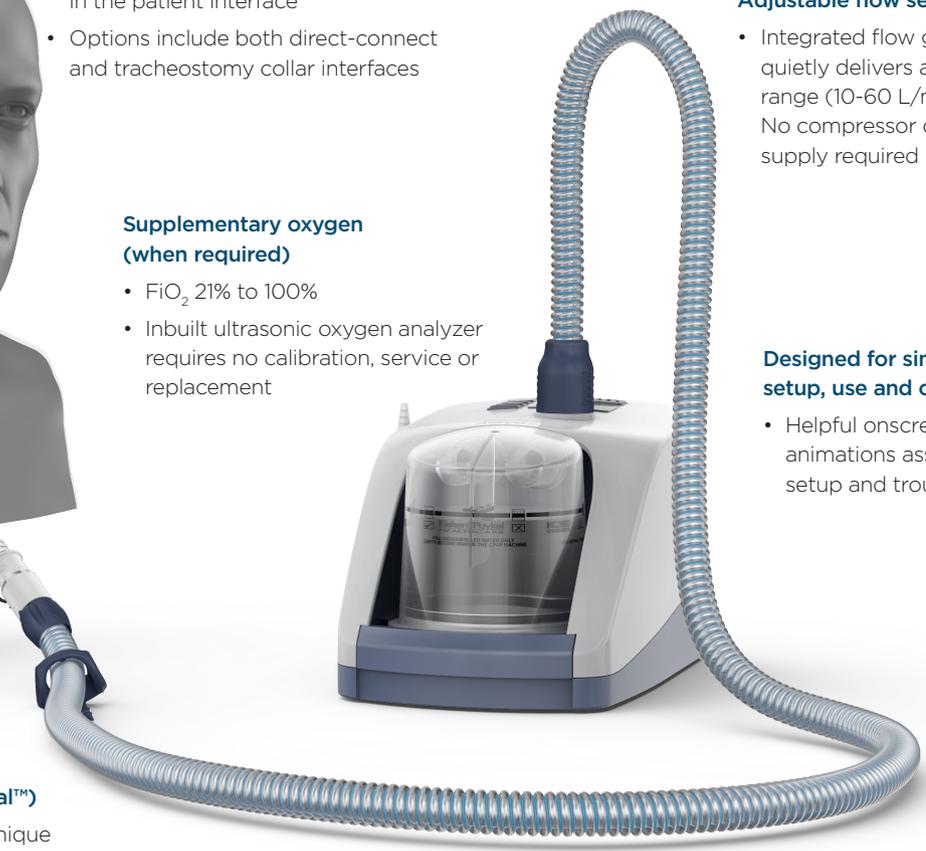

Water vapor vs. droplets

The F&P myAIRVO™ 2 humidifier is designed to generate water vapor (molecules).

- Nebulizers and other aerosol delivery systems are designed to generate water droplets. Unlike water vapour, the droplets are large enough to carry pathogens.¹

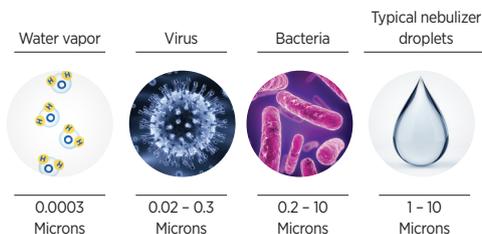

Product codes	
PT100xx	myAIRVO 2 humidifier
900PT560	AirSpiral Heated breathing tube (10-pack)
900PT560E	AirSpiral Heated breathing tube (1-pack)
OPT970	Optiflow+ Tracheostomy interface (20-pack)
OPT970E	Optiflow+ Tracheostomy interface (1-pack)
OPT971	Sputum guard spare part (20-pack)
OPT980	Mask interface adapter (20-pack)
OPT980E	Mask interface adapter (1-pack)
900PT290E	MR290 auto-fill water chamber
HC360	HC360 reusable water chamber
900PT400	Compact stand
900PT401	Water bag (2-pack)

1. Takigawa K, Fujita J, Negayama k, Yamagishi Y, Yamaji Y, Ouchi K. Nosocomial Outbreak of Pseudomonas cepacia Respiratory Infection in Immunocompromised Patients Associated with Contaminated Nebulizer Devices. Jpn J Infect Dis. 1993; 67(11):1115-1125.

Optiflow, Evaqua, AirSpiral, myAIRVO, are trademarks of Fisher & Paykel Healthcare Limited. For patent information, see www.fphcare.com/ip

Find out more at www.fphcare.com/myairvo

Fisher & Paykel
HEALTHCARE