

High-performance seal
enables Bridge-free NIV

A design that delivers on comfort,
stability and performance.



Dynamic Support Technology™

The F&P Visairo™ mask features the Floating Seal and Stability Wings which work in synergy to create Dynamic Support Technology - a design that delivers on comfort, stability and performance.




96% of caregivers found the mask performance to be good or very good.¹



Stability Wings

The Stability Wings in the frame provide support for the Floating Seal to remain comfortable and stable on your patient's face.

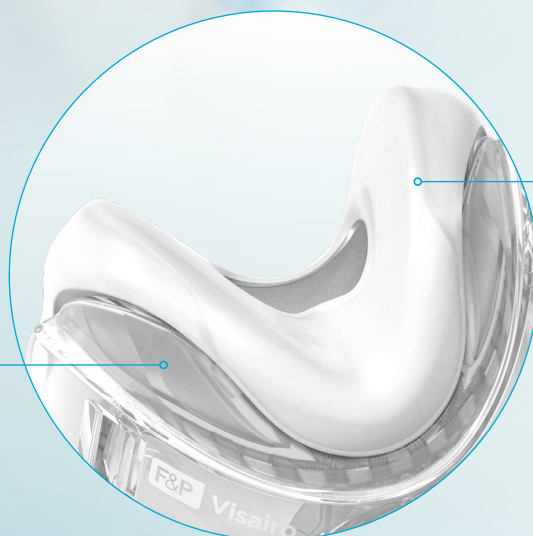


Floating Seal

The soft, silicone Floating Seal is decoupled from the mask frame so that it floats comfortably under your patient's nose.


Bridge-free NIV experience

Dynamic Support Technology
delivers a Bridge-free NIV experience.

The F&P Visairo mask removes pressure
from the bridge of the nose and clears the
way for an open field of vision, without
compromising on stability and performance.

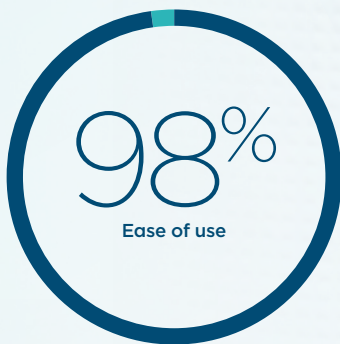

98% of caregivers found
the Visairo mask to be easy
or very easy to use.³

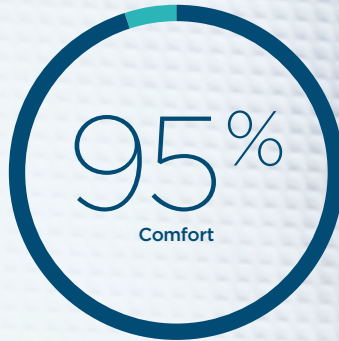

95% of caregivers found
perceived patient comfort
to be good or very good.²





Three mask sizes

The F&P Visairo mask is available in three sizes: A, B and C, which correspond to the dimensions of your patient’s nose.



TubeFit™ seal zones

The seal zones, featured on both sides of the mask, are designed to seal around nasogastric and orogastric tubes. This helps reduce the need to over-tighten the headgear.



Three mask models

The The F&P Visairo mask is available in non-vented and vented options for single or dual limb set-ups.

F&P Visairo Full Face Bridge-free Mask

Mask Version	Product Order Code
Non-vented mask for single-limb circuit	RT075A
	RT075B
	RT075C
Non-vented mask for dual-limb circuit	RT076A
	RT076B
	RT076C
Vented mask for single-limb circuit	RT077A
	RT077B
	RT077C

1. 89 out of 93 responses from caregivers rated the Visairo mask performance as good or very good. Internal validation trials conducted with 100 feedback forms filled out by caregivers from 12 different hospitals across the United States of America, Canada, and Australia.

2. 90 out of 95 responses from caregivers rated the perceived patient comfort on the Visairo mask as good or very good. Internal validation trials conducted with 100 feedback forms filled out by caregivers from 12 different hospitals across the United States of America, Canada, and Australia.

3. 93 out of 95 responses from caregivers rated the Visairo mask as easy or very easy to use. Internal validation trials conducted with 100 feedback forms filled out by caregivers from 12 different hospitals across the United States of America, Canada, and Australia.

F&P, Bridge-free, Dynamic Support Technology, TubeFit, and Visairo are trademarks of Fisher & Paykel Healthcare.
For patent information, please see www.fphcare.com/ip. PM-622647 REV D © 2023 Fisher & Paykel Healthcare Limited